Silk'n™
Beautiful Technology

Cleansing Box for Glide and Infinity

INSTRUCTIONS FOR USE

EN SV

DE DA

FR BG

NL HU

ES RO

PT CS

GR SL

TR PL

IT RU

NO

Signs and Symbols

	Products marked with this symbol fulfil the requirements of the EU directives.
IP30	The device is protected against the ingress of objects larger than 2.5mm diameter, such as fingers, to protect against electric shock.
	Do not use in wet environments.
	Do not dispose with household waste. Dispose in accordance with your country's laws and help protect the environment. In accordance with the WEEE directive.

User Manual

This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.

Children should not play with the appliance.

Cleaning and user maintenance should not be undertaken by children without supervision.

Dispose of the packaging at a recycling point which sorts materials by type.

Products marked with this symbol should not be disposed of along with household waste! You are legally bound to dispose of old appliances separately from household waste.


Do not dispose of batteries with household waste. You are legally bound to dispose of all batteries at a local authority collection point or to return them to a battery retailer.

Information about where old appliances can be disposed of free of charge is available from your local authority.

Silk'n Cleansing Box device description

1. Disinfection heads
2. Device activated on/off button
3. Internal blue LEDs
4. External blue LEDs
5. Release lip
6. Battery compartment

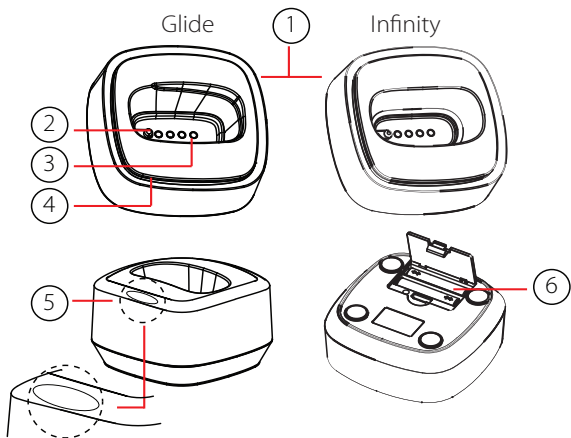Intended use

Devices which work in direct contact with the skin need to be thoroughly disinfected after each treatment. Silk'n Cleansing Box is a small, lightweight and easy to use at-home disinfectant which will help keep your device clean and safe to use at all times!

Package contents

Cleansing Box device
2x1.5V LR03 (AAA) batteries

Instructions for use

Batteries

- Use batteries with the correct polarity. Do not mix new and used batteries.
- Remove exhausted batteries from the device. If the appliance is to be stored unused for a long period, please remove the batteries.
- The supply terminals are not to be short circuited.
- Do not recharge batteries.

How to use

Refer to the diagram to see the relevant parts.

1. Place the batteries inside the compartment (#6).
2. Select the correct disinfection head for your device (#1). To swap heads, release the head from the device. Use a small blunt tool to lever the head off the device carefully. There is a small lip (#5) to insert the tool. Make sure that the correct head is placed in the direction shown in the diagram.
3. Next, place the device into the disinfection head, with the treatment surface face down. It will fit snugly if the correct head is in place.
4. Push the device down so it presses on the On/Off button to activate the box. The blue LEDs will turn on automatically for 5-6 minutes. This is the amount of time required for a full disinfection process. It is also possible to stop the disinfection process before time, by pushing the device down again so it presses the On/Off button and deactivates the box.

Warning! We recommend that you do not operate the On/Off switch manually to protect your eyes from concentrated light.

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Customer Service

See our website for more information: www.silkn.eu. This manual is also available as a PDF download from the website. Contact Customer Service to report unexpected operations or events, or any other problem with your device. Our contact details are:

Country	Service number	Email
Belgium	0900-25006	serviceconsommateurbe@silkn.eu
France	0891-655557	serviceconsommateurfr@silkn.eu
Germany	089 51 23 44 23	kundenservice@inno-essentials.de
Italy		servizioclientiit@silkn.eu
Netherlands	0900-2502217	klantenservicentl@silkn.eu
Spain	807088200	servicioalcliente@silkn.eu
UK	0906-2130009	customerareuk@silkn.eu
Other countries	+31 (0)180-330 550	info@silkn.com

PB111475A

A stylized, grey, wavy graphic element resembling a drop or a wave, positioned above the brand name.
HomeSkinovations™
info@silkn.com
www.silkn.eu